## InnovaGoods®

### User's Manual

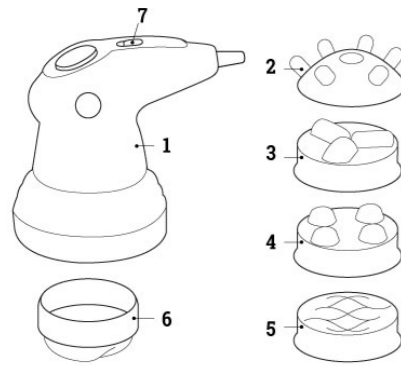
Fig. 1


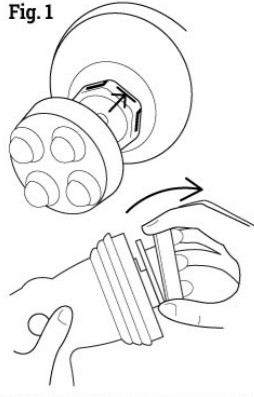
Fig. 2
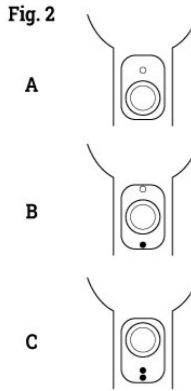
## ENGLISH

### TECHNICAL SPECIFICATIONS

- Voltage: 220-240 V
- Frequency: 50 Hz
- Power: 25 W
- Speed control: up to 2500 r.p.m.

- Infrared light

### COMPONENTS

1. Suction & Heat Anti-Cellulite MassagerInnovaGoods®
2. Finger pressure head (legs and glutes): deep massage, calming, cellulite
3. Roller head (legs and glutes): draining massage, cellulite, flaccidity
4. Ball head (abdomen and arms): firm skin and muscles
5. Undulating head (whole body): relaxing and toning massage
6. Protective fabric mesh: for undulating head in areas with hair
7. Infrared heat: thermo-active effect during the massage to mobilise the accumulated fat and toxins

### INSTRUCTIONS

- Connect the Suction & Heat Anti-Cellulite Massager InnovaGoods® to the mains electricity.
- Insert the desired head (figure 1) and move the switch (figure 2): Position A: to switch the massager on /Position B: to activate the infrared heat / Position C: to switch off.
- Rotate the potentiometer to adjust the speed.
- Apply essential oils and moisturising cream to the skin which is to be treated so it helps to minimise the skin's redness (never when using the protective mesh).
- Hold the massager with both hands and make slow and circular movements on the area to be treated; start with a slow speed and increase according to your needs. Do not press too hard, just slide gently.
- Recommended daily use: between 15-20 mins (5 mins in each area).
- Clean the heads after each use.
- You can apply an anti-cellulite or firming cream after using the device to boost results.

### WARNINGS

Do not use the infrared heat function on varicose veins. Do not use the massager on the abdomen if you are pregnant. NEVER massage on a bone (knee, ankle or elbow). Make sure that the voltage of the power source corresponds to that of the appliance. Unplug the appliance whenever you are not using it. Do not pull on the cable and ensure it does not tangle. Keep the cable away from hot and sharp surfaces. Before cleaning, make sure that appliance is disconnected from the power supply and has cooled completely. Clean with a damp cloth. Do not use abrasive products for cleaning. Do not immerse the appliance or cable in liquid. Do not handle it with wet hands Do not use near water or damp areas. Suitable for domestic use only. Use the device for its intended purpose only. Store the device in a dry, dust-free place. Do not leave the appliance unsupervised when it is turned on. This appliance is not a toy. Keep out of reach of children. Do not modify the appliance or any of its components. If the appliance, the cable or the pieces looked damaged or do not work correctly, do not use them or attempt to repair them; they should only be repaired by a qualified technician. This device is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by an adult responsible for their safety. This device is not household waste. Please follow the environmental

regulations in place for its recycling Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

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